

$\alpha\beta$ Axis Extended Contact Ball Slide Goniometers


Stage Size 40 × 40 mm

GOHT-40B

RoHS

High precision goniometers with integrated bearing ways for superior stiffness, accuracy and durability.

- Centering accuracy 10x better than typical dovetail goniometers.
- Designed for use in a variety of orientations. Not limited to horizontal mounting.
- Integral-guide structure offers high performance of rotation center displacement $\phi 0.007\text{mm}$.
- Disc clamp improves the stability of clamping.
- Designed for use with precision alignment pins for ease of assembly.

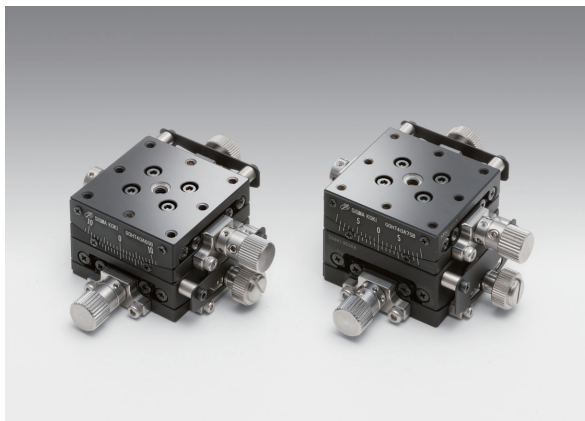


Guide

- Attachment pins (GOHT-AP) can be purchased separately.
[Reference](#) E010 [Catalog Code](#) W7112
- Contact our Sales Division for information about custom center heights.
- Motorized version (OSMS-40B***) also available.
[Reference](#) G094 [Catalog Code](#) W9098



OSMS-40B***

- For the surface finishing, the white coating of electroless nickel plating is also available. [WEB Reference](#) [Catalog Code](#) W7104

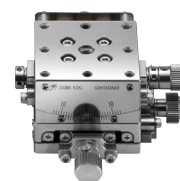

Attention

- Locking mechanism is operated using thumbscrew on side of stage.

Specifications

Part Number	METRIC	GOHT-40B20B	GOHT-40B40B	GOHT-40B60B
	INCH	GOHT-40B20BUU	GOHT-40B40BUU	GOHT-40B60BUU
(Opposite Model)	METRIC	GOHT-40B20BR	GOHT-40B40BR	GOHT-40B60BR
	INCH	GOHT-40B20BRUU	GOHT-40B40BRUU	GOHT-40B60BRUU
Stage Size [mm]		40×40	40×40	40×40
Axes of Travel		$\alpha\beta$ axis		
Rotation Center Height [mm]		20±0.2	40±0.2	60±0.2
Angle Range [°]	(Upper) β axis	±20	±10	±8
	(Lower) α axis	±10	±8	±6
Vernier Readable Resolution		0.1°	0.1°	0.1°
Travel per Knob Rotation [°]	(Upper) β axis	about 2.34	about 1.41	about 1.08
	(Lower) α axis	about 1.41	about 1.08	about 0.89
Rotation Center Displacement [mm]		$\phi 0.007$ or less	$\phi 0.007$ or less	$\phi 0.007$ or less
Guide Method		Extended contact ball bearing guide		
Primary Material		SUS440C quench hardened		
Finish		Super black chrome		
Load Capacity [N]		58.8 (6.0kgf)	58.8 (6.0kgf)	58.8 (6.0kgf)
Max. Moment Capacity	Pitch [N·m]	2.0	2.0	2.0
	Roll [N·m]	2.0	2.0	2.0
	Yaw [N·m]	2.0	2.0	2.0
Parallelism [μm]		30	30	30
Weight [kg]		0.4	0.4	0.4

Application Systems

Optics & Optical Coatings

Opto-Mechanics

Bases

Manual Stages

Actuators & Adjusters

Motorized Stages

Light Sources & Laser Safety

Index

Guide

X Axis Stages

XY Axis Stages

Z Axis Stages

XZ Axis Stages

XYZ Axis Stages

Rotation Stages

Goniometer

Tilt Stages

Vacuum

Ball Bearing Guide

Crossed Roller

Dovetail

Lapping

V Groove Screw

Others

15 × 15 mm

25 × 25 mm

40 × 40 mm

60 × 60 mm

65 × 65 mm

80 × 80 mm

100 × 100 mm

120 × 120 mm

Others

