

HIT Controller Package

HIT-PKG

RoHS **Catalog Code** W9048

This product combines the HIT series extensible stage controllers into a single package.

- This product combines one HIT-MV (Extensible Stage Controller (Master)) and two to a maximum of eight extensible stage controller slaves (HIT-SA) into a single package.
- The power supply for each HIT has been unified (AC100V-120V/200-240V), and the link cables are internally routed, tidying up the usually complex cable arrangement.
- All axes can be carried together at once.



HIT-PKG-08
FRONT

BACK

Common Specifications

Power Voltage	AC100~240±10%[V] 50/60[Hz]
Power Consumption	750VA or lower
Dimensions (W×D×H)	450×304×364mm (excluding protrusions)
Storage Temperature	-20 - +60°C
Operating Temperature	+5 - +35°C
Ambient Humidity	20 - 80%RH (without condensation)
Primary material	Aluminum
Finish	Black semi-gloss paint

Specifications

Part Number	HIT-PKG-02	HIT-PKG-03	HIT-PKG-04	HIT-PKG-05	HIT-PKG-06	HIT-PKG-07	HIT-PKG-08
JP Yen [¥]	340,000	400,000	460,000	610,000	670,000	730,000	790,000
Number of Control Axes	2	3	4	5	6	7	8
Package Contents							
HIT-MV	1	1	1	1	1	1	1
HIT-SA	2	3	4	5	6	7	8
Weight [kg]	9.0	9.5	10.0	16.5	17.0	17.5	18.0



Guide

- Sample programs are available for download on our website.
 - SG Sample 32/64 bit version for Windows® (RS232C (Virtual COM), for TCP/IP)
 - LabVIEW for RS232C (Virtual COM), for TCP/IP)
- If you wish to include the HIT-SA-M2 (slave for 2-phase stepping motor) in your configuration, please contact our sales department.

Attention


- Please use the included power cable.
- As this is a custom-made product, modifications such as adding more axes cannot be accommodated.

Outline Drawing


(Alt : mm)

HIT-PKG-02 to 04

HIT-PKG-05 to 08

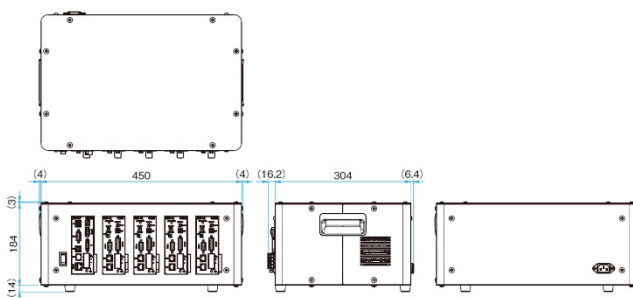

The package contents vary depending on the controlled axis (common for HIT-PKG-02 to 04).

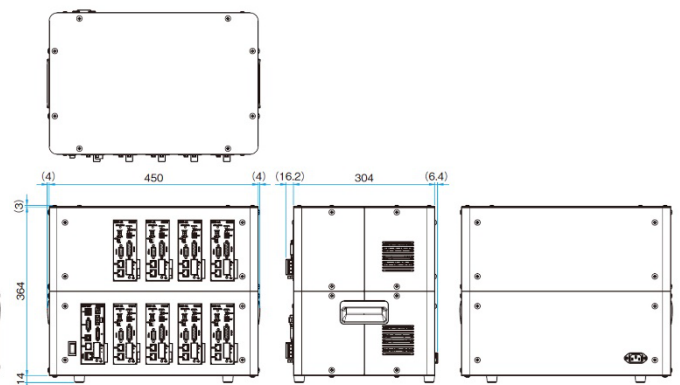

The package contents vary depending on the controlled axis (common for HIT-PKG-05 to 08).



Tokyo Head Office
5F, SIGMA KOKI Tokyo Head office 1-19-9,
Midori, Sumida-ku Tokyo 130-0021 JAPAN

Osaka Branch
4-9-28 Nishi-Nakajima,
Yodogawa-ku, Osaka 532-0011

Kyushu Sales Office
1-17-25 Hakata-eki-higashi, Hakata-ku,
Fukuoka 812-0013

INTERLAB Sp. z o.o. | Instrumenty badawcze i pomiarowe

Kosiarzy 37 paw. 20, 02-953 Warszawa | +48 22 840 81 80 | interlab@interlab.pl | interlab.pl

Extensible Stage Controller | HIT-MV/SA

RoHS

CE

Modular stage control system with Master controller and from one to eight slave Axes.

- Control with RS232C/USB/Ethernet interfaces is available.
- Full closed loop control is possible when used in positioning stage equipped with encoder.



Guide

- ▶ Sample programs are available for download on our website.
- SG Sample 32/64 bit version for Windows® (RS232C (Virtual COM), for TCP/IP)
- LabVIEW for RS232C (Virtual COM, for TCP/IP)

Attention

- ▶ The power specifications for HIT-MV are DC +24V 1A. The power specifications for HIS-SA are DC +24V 2A. Depending on the number of axes in the stage, please prepare a power supply ranging from 3A (for 1 axis) to 17 axes (for 8 axes). Please purchase the AC-ADP-2427 (AC adapter) or prepare a power supply separately.

Part Name	Part Number
Extensible Stage Controller (Master)	HIT-MV
Extensible Stage Controller (Slave)	HIT-SA
LAN Cable	LAN-2
AC Adapter	AC-ADP-2427

Primary Functions

Part Number	HIT-MV	HIT-SA
Controller Function	○	—
Number of Control Axes	8	—
Stored Program Control	○	—
Feedback Control	—	OSMS (CS) series HST (GS) series
Circular Interpolation Control	○	—
Linear Interpolation Control	3 axes	—
Driver Function	—	Micro-step
Micro-step (Max. Division)	—	250
Driving Current (A/phase)	—	0.25 – 1.4

General Specifications

Part Number	HIT-MV	HIT-SA
Power Voltage	DC24V 1A	DC24V 2A
Power Consumption	24VA	48VA
Operating Temperature	5 – 40°C	
Storage Temperature	–20 – 60°C	
Ambient Humidity	20 – 80%RH (without condensation)	
External Dimensions (W×H×Dmm)	130×120×50	130×120×50
Weight (kg)	0.62	0.62

Interface

GP-IB	—
RS232C	○
USB	○*
Ethernet	○

*Acknowledge as virtual serial port from external PC.

Performance Specifications

Coordinate Indication Range	—
Max. Travel to Set	134,217,727
Max. Driving Speed (pps)	500,000
Min. Driving Speed (pps)	1
Acceleration/Deceleration Time (ms)	1 – 1,000

I/O Specification

Origin Sensor	○
Proximity Sensor	○
CW (+) Limit	○
CCW (–) Limit	○
General Purpose Input	4 points
General Purpose Output	4 points
Control Input	—
Control Output	—
Trigger Output	—

Control Command

Machine Origin Return	○
Theoretical Origin Setting	○
Relative Position Drive	○
Absolute Position Drive	○
Jog Operation	○
Position Appointment	—
Circular Interpolation Control	○
Linear Interpolation Control	○
Drive	○
Deceleration Stop	○
Emergency Stop	○
Speed Setting	○
Motor Free/Hold	○
Port Input	○
Port Output	○

Optional

CJ-200A	—
JS-301	—
JB-401	—
JD-101	—
SJT-02	—
MD-400	—

Application Systems

Optics & Optical Coatings

Opto-Mechanics

Bases

Manual Stages

Actuators & Adjusters

Motorized Stages

Light Sources & Laser Safety

Index

Guide

Controllers/Drivers

Softwares

Stepping Motor

AC Servo Motor

Cables

Piezo

X Translation

Theta Rotation

Goniometer

Vacuum

Options

40 × 40 mm

60 × 60 mm

80 × 80 mm


85 × 85 mm

100 × 100 mm

120 × 120 mm

Others

■ HIT System Chart

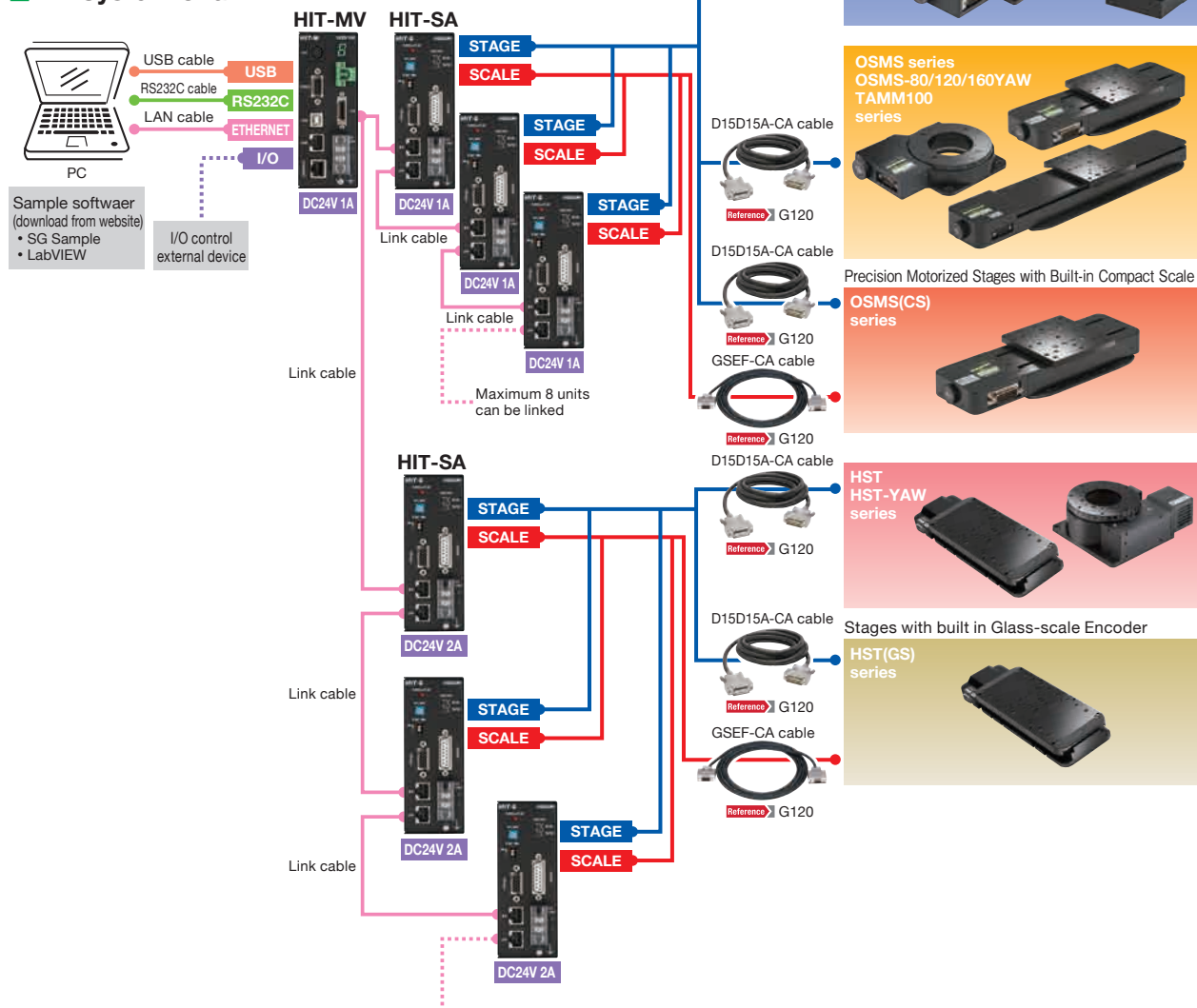

Stepping Motor Stage

OSMS-A/B
OSMS-ZF
OSMS-40/60YAW
SGSP-ACT-B0
TAMM40/60
HPS
HDS
series



OSMS series
OSMS-80/120/160YAW
TAMM100
series



Precision Motorized Stages with Built-in Compact Scale

OSMS(CS) series




HST
HST-YAW
series



Stages with built in Glass-scale Encoder

HST(GS)
series

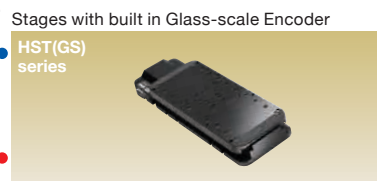

Application Systems

Optics & Optical Coatings

Opto-Mechanics

Bases

Manual Stages

Actuators & Adjusters

MotORIZED STAGES

Light Sources & Laser Safety

Index

Guide

Controllers/Drivers

Softwares

Stepping Motor

AC Servo Motor

Cables

Piezo

X Translation

Theta Rotation

Goniometer

Vacuum

Options

11

40 × 40 mm

60 × 60 mm

80 × 80 mm

132 132

120 × 120 mm

Others